


# INSTRUMENT DESIGN & MFG. CO.

## HYDRO ASPIRATOR CLEANING INSTRUCTIONS

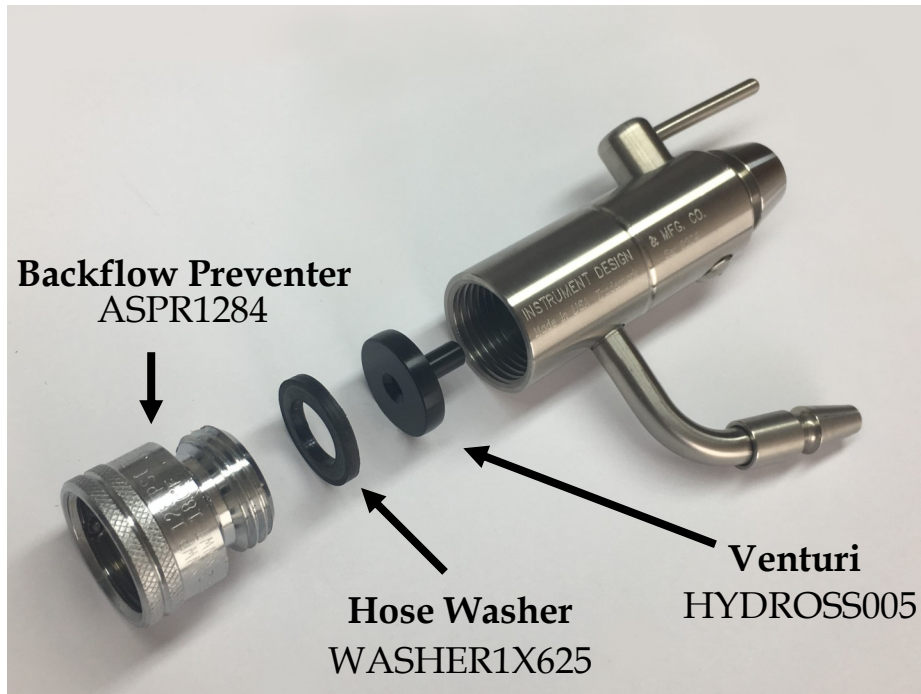

[Click here for a video on disassembling the Vacua Hydro Aspirator™](#)

### **DO NOT**

- Stick implements into venturi hole
- Disassemble backflow preventer
- Remove aspirator valve (unless it needs repair)

### Disassembly

1. Remove the **backflow preventer** by unscrewing it from the aspirator.
2. Pry the **hose washer** out of the aspirator body.
3. Tap aspirator bottom to remove **venturi**.
4. Run water through aspirator to remove any debris.

### Reassembly

1. Slip **venturi** back into aspirator. If it doesn't go straight in, try again.
2. Press the **hose washer** down against the **venturi**.
3. Screw **backflow preventer** back onto aspirator.

Supplemental Products: Hydro Aspirator Repair Kit (HYDROREKIT) includes venturi, hose washer, 2 o-rings and grease (for valve).