


INSTRUMENT DESIGN & MFG. CO.

VACUA HYDRO ASPIRATOR™ INSTRUCTION SHEET

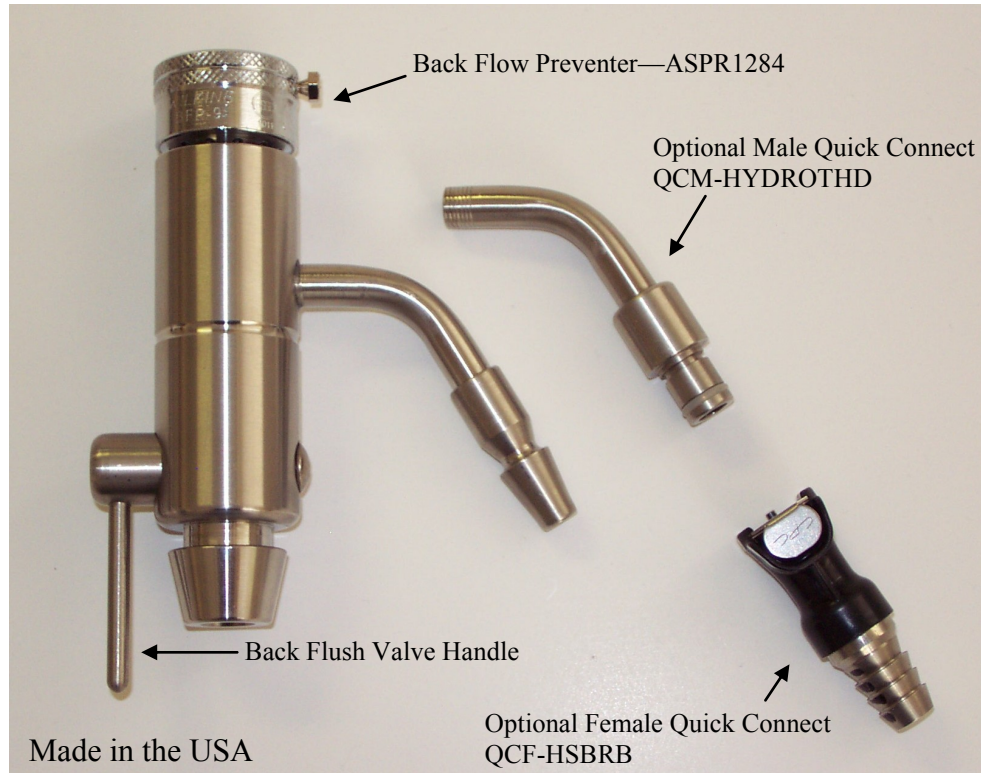

INSTALLATION: To install simply attach to any sink faucet with garden hose type threads. Attach suction hose to 90 degree elbow on side of unit. Next attach a splash tube to outlet at the base of the Hydro Aspirator.

OPERATION: To operate place the back flush valve handle in the vertical position and turn on the water. Vacuum should be present at end of inlet hose. Increase water flow to increase vacuum. The Hydro Aspirator does not require water pressure and may not function properly if faucet is not turned on sufficiently. Adjust water flow to obtain the level of vacuum required. Excessive vacuum may draw excessive amounts of material such as clots in to the Hydro Aspirator potentially clogging the Hydro Aspirator. If this occurs try reducing the water flow or using a Trocar with smaller holes.

CLEANING: To clear a blockage in your Hydro Aspirator you can turn the back flush handle to the horizontal position and force water backward through the inlet hose (Note: Some water will still come out the base of the Hydro Aspirator, this is normal). In the event material becomes clogged and disassembly is required, simply unscrew the Hydro Aspirator from the Back flow preventer and remove the rubber washer inside. Tap unit upside down to remove the venturi (Black cone shaped piece inside DO NOT Stick anything into the venturi, Only water passes through it and putting anything in it will damage it . A replacement Venturi is available from your distributor. #HYDROSS005). Clean any obstruction from the inside of the Hydro Aspirator and reassemble.

INSTRUMENT DESIGN & MFG. CO.

25 Harvey Road
Unit 6
Bedford, NH 03110

Phone: 603-606-5412
Fax: 603-606-5101
E-mail: injector@inj.mv.com
Web Site: www.embalminginstrument.com